

SHZ-95B circulating water vacuum pump

High temperature resistant, anti-corrosion and good insulating property
Stepless speed regulation, stable in low speed and powerful in high speed

SHZ-95B循环水式真空泵

耐高温，耐腐蚀，绝缘性能好
无极调速，低速平稳，高速强劲



☎ whatsapp: +86-13783661974

☎ : +86-371-67447888

☎ : +86-371-67447999

🔍 : www.rotovap.cn

@ : rotovap@lanphan.com

PRODUCT DESCRIPTION

产品描述

SHZ-95B Vertical Circulating Water Vacuum Pump is a new type of vacuum pump designed with circulating water as working fluid and negative pressure generated by jet principle. It is not only a vacuum pumping device, but also provides the reaction device. Circulate cooling water. The pump has the characteristics of no oil, no pollution, corrosion resistance, low noise, convenient and flexible. It adopts anti-corrosion movement or stainless steel movement, vertical motor and multi-purpose machine.

SHZ-95B Vertical Circulating Water Vacuum Pump Features

- (1) This machine adopts five taps, which can be used separately or in parallel, and is equipped with a vacuum gauge;
- (2) Non-polluting, corrosion-resistant, low noise, convenient and durable, saving water and electricity;
- (3) The main machine is made of two models: stainless steel movement and anti-corrosion material movement;
- (4) Small size, light weight, convenient movement and beautiful appearance.
- (5) The machine adds water level and the indicator light is displayed; high, medium and low.

SHZ-95B 型立式循环水真空泵以循环水为工作流体，采用射流原理产生负压而设计的一种新型真空抽气泵，它不仅是一种真空抽气装置，同时还能向反应装置中提供循环冷却水。该泵具有不用油、无污染、耐腐蚀、噪音低、方便灵活等特点。采用防腐机芯或不锈钢机芯，立式电机，一机多用。

SHZ-95B 型立式循环水真空泵产品特点

- (1) 本机采用五抽头，可单独或并联使用，装有一个真空表；
- (2) 无污染、耐腐蚀、噪音低、方便耐用、省水省电；
- (3) 主机采用不锈钢机芯和防腐材质机芯两种型号制造；
- (4) 体积小，重量轻，移动方便，外型美观。
- (5) 本机加水位，指示灯显示；高、中、低。

Product parameters

产品参数

■ 循环水立式真空泵 /Circulating water vertical vacuum pump

型号 Model	SHZ-95B
流量 Flow L/min	80L/min
扬程 Lift	12M
功率 Power	370W
单头抽气量 Single head pumping	10L/min
储水箱容积 Tank volume	50L
真空度 Vacuum degre	-0.098Mpa
电源 Power Supply	220V/110V 60Hz/50Hz
抽头个数 Draw a number	Single meter 5 taps单表5抽头

Product Details

产品详情




Wear-Resisting Rubber Wheel

Equipped with wear and pressure resistant rubber wheel, which is convenient to move and long service life

Vertical Circulating Water Vacuum Pump

Five taps of single or multiple using. Parallel connected with five-way pipe (optional) has larger sucking rate, beable to satisfy needs of large scale rotary evaporator or reactors

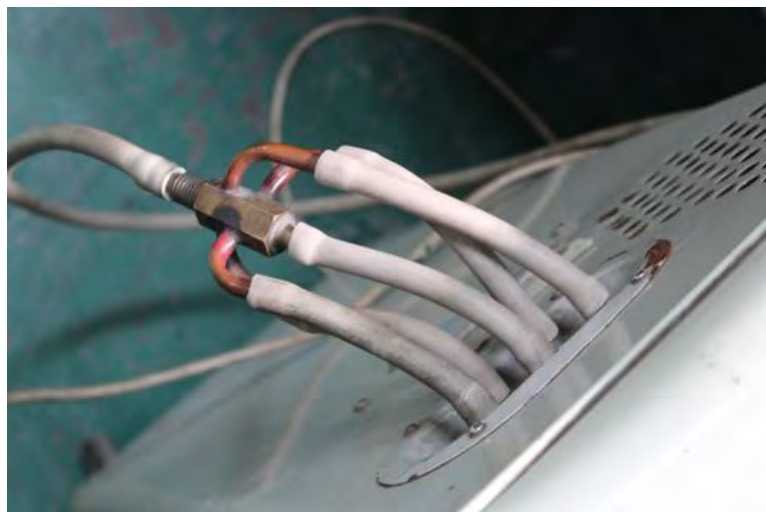

Inlet/Outlet

Clear label on water inlet and outlet, details manifesting excellent services.




Product Details

产品详情

Vacuum Meter

Vacuum meter makes experiment conditions clear at a glance.

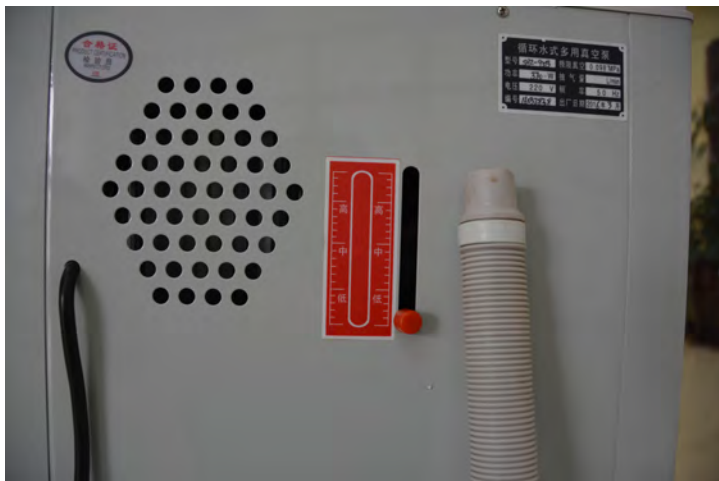

Vertical Circulating Vacuum Pump Water Level Display

High, middle and low water level display, convenient and easy to control experiment condition.

Internal Structure

Considered form inside and outside reflecting manner of big factory.



MATCHING PRODUCTS

配套产品



SHZ-95B circulating
water vacuum pump
SHZ-95B循环水式真空泵




Rotary evaporator
旋转蒸发器



Reactor
反应釜



QUALITY CERTIFICATION

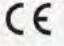合格证书

 ISET S.r.l. <small>Sede Legale e Uffici: Via Donato di Sanguis 51 40024 Moglia (MO) Tel. e fax +39 (0)76 589863 jset@set-italia.com www.set-italia.com</small> <small>Cap. Soc. i.v. Cod. Fisc. e P.IVA Reg. Imprese REA</small>		<small>€ 10.200,00 02 232 750 369 02 232 750 369 MN 0221088</small>
CERTIFICATE OF COMPLIANCE Certificado de Conformidade - Сертификат соответствия - Konformitätserklärung		
D) APPLICANT: Henan Lanphan Trade Co., Ltd. Room 801, Building B, North Jianshe Road, West Qingling Haod Road, Zhongyuan District, Zhengzhou City, China MANUFACTURER: Henan Lanphan Industry Co., Ltd. North Jianshe Road, West Qingling Road, Zhongyuan District, Zhengzhou City, China	3) CERTIFICATE NO.: IT0524LP16011801 TEST REPORT(S) NO.: DL-2018010006E DL-2018010009S DL Certification & Testing Co., Ltd.	WITH REFERENCE TO EC DIRECTIVE APPLIED: Low Voltage Directive 2014/35/EU Electromagnetic Compatibility 2014/30/EU HARMONIZED STANDARDS APPLIED: EN 61326-1:2013, EN 61000-3-2:2014, EN 61000-3-3:2013, EN 61010-1:2010
5) PRODUCT CHARACTERISTICS: Circulating Water Vacuum Pump TRADE MARK:  MODEL(S): SHZ-95B, SHZ-D111, SHZ-GA, SHZ-CB, SHZ-C, SHZ-CD, SHZ-2000		
REMARK: The verification has been carried out on a voluntary basis. We attest that a TCF is in place. The product(s) satisfies the requirements of the Certification Mark of ISET, in reference to the above list standard(s). The above compliance mark can be fixed on the product(s) according to the ISET regulation about its release. This verification doesn't imply assessment of the production and the products. Notice of the CE marking: The label of the CE marking: Not less than 5mm height. Before putting the product(s) into market, CE marking and EC declaration are duties of the manufacturer. The manufacturer is responsible to start the CE marking certification procedure and to perform the activities according to all the relative directive(s). 		
Q) DATE OF ISSUE: 16/01/2018 DATE OF EXPIRE: 15/01/2021 CERTIFICATION MANAGER:  		

 ISET S.r.l. <small>Sede Legale e Uffici: Via Donato di Sanguis 51 40024 Moglia (MO) Tel. e fax +39 (0)76 589863 jset@set-italia.com www.set-italia.com</small> <small>Cap. Soc. i.v. Cod. Fisc. e P.IVA Reg. Imprese REA</small>		<small>€ 10.200,00 02 232 750 369 02 232 750 369 MN 0221088</small>
CERTIFICATE OF COMPLIANCE Certificado de Conformidade - Сертификат соответствия - Konformitätserklärung		
D) APPLICANT: Henan Lanphan Trade Co., Ltd. Room 801, Building B, North Jianshe Road, West Qingling Haod Road, Zhongyuan District, Zhengzhou City, China MANUFACTURER: Henan Lanphan Industry Co., Ltd. North Jianshe Road, West Qingling Road, Zhongyuan District, Zhengzhou City, China	3) CERTIFICATE NO.: IT0524LP16011801 TEST REPORT(S) NO.: DL-2018010006E DL-2018010007S DL Certification & Testing Co., Ltd.	WITH REFERENCE TO EC DIRECTIVE APPLIED: Low Voltage Directive 2014/35/EU Electromagnetic Compatibility 2014/30/EU HARMONIZED STANDARDS APPLIED: EN 61326-1:2013, EN 61000-3-2:2014, EN 61000-3-3:2013, EN 61010-1:2010
5) PRODUCT CHARACTERISTICS: Rotary Evaporator TRADE MARK:  MODEL(S): RE-501, RE-201D, RE-52A, RE-52C, RE-2010, RE-2000A, RE-2000B, RE-2000E, RE-301, RE-3000A, RE-3000B, RE-3000E, RE-500X, RE-1002, RE-2001Z, RE-3002, RE-5003, RE-100X, R-1003, R-1010, R-1020, R-1050		
REMARK: The verification has been carried out on a voluntary basis. We attest that a TCF is in place. The product(s) satisfies the requirements of the Certification Mark of ISET, in reference to the above list standard(s). The above compliance mark can be fixed on the product(s) according to the ISET regulation about its release. This verification doesn't imply assessment of the production and the products. Notice of the CE marking: The label of the CE marking: Not less than 5mm height. Before putting the product(s) into market, CE marking and EC declaration are duties of the manufacturer. The manufacturer is responsible to start the CE marking certification procedure and to perform the activities according to all the relative directive(s). 		
Q) DATE OF ISSUE: 16/01/2018 DATE OF EXPIRE: 15/01/2021 CERTIFICATION MANAGER:  		

 CERTIFICATION EUROPE HONG KONG	
Certificate of Registration Certificate NO: 10134302E This is to Certify that the Environmental Management System of Henan Lanphan Industry Co., Ltd. Room 801, Block B, North of Jianshe Road, West of Qingling Road, Zhongyuan District, Zhengzhou City, He'nan Province, China Has been audited to the following Environmental Management System standard: ISO14001:2015 This system is valid for the Sales of Valves, pipe fittings, rubber products, plastic products, metal products, glassware, water treatment equipment, stainless steel products, children's toys, needles textiles, mechanical and electrical products, amusement equipment, building materials, machinery equipment, hardware tools; Design and construction of anti-corrosion insulation engineering, municipal engineering, highway engineering; Construction labor subcontracting; Electromechanical equipment installation engineering; Installation of water and electricity; Integrated wiring system engineering; Installation and maintenance of computer; Technical development, technology transfer and technical consultation of computer software; Technical development of photoelectric integration; Property services; Cleaning services involved in environmental management activities Date of issue: Jan. 09, 2018 Date of expiry: Jan. 08, 2021 This certificate will not remain valid only if the certified organization accepts at least one surveillance audit annually within the validity period of the certificate in which the surveillance audit conforming mark is in the designated position on the certificate.    Issued by:  CERTIFICATION EUROPE(HONGKONG)LIMITED WWW.CCE-HK.COM	

 CERTIFICATION EUROPE HONG KONG	
Certificate of Registration Certificate NO: 10134302Q This is to Certify that the Quality Management System of Henan Lanphan Industry Co., Ltd. Room 801, Block B, North of Jianshe Road, West of Qingling Road, Zhongyuan District, Zhengzhou City, He'nan Province, China Has been audited to the following Quality Management System standard: ISO9001:2015 This system is valid for the Sales of Valves, pipe fittings, rubber products, plastic products, metal products, glassware, water treatment equipment, stainless steel products, children's toys, needles textiles, mechanical and electrical products, amusement equipment, building materials, machinery equipment, hardware tools; Design and construction of anti-corrosion insulation engineering, municipal engineering, highway engineering; Construction labor subcontracting; Electromechanical equipment installation engineering; Installation of water and electricity; Integrated wiring system engineering; Installation and maintenance of computer; Technical development, technology transfer and technical consultation of computer software; Technical development of photoelectric integration; Property services; Cleaning services Date of issue: Jan. 09, 2018 Date of expiry: Jan. 08, 2021 This certificate will not remain valid only if the certified organization accepts at least one surveillance audit annually within the validity period of the certificate in which the surveillance audit conforming mark is in the designated position on the certificate.    Issued by:  CERTIFICATION EUROPE(HONGKONG)LIMITED WWW.CCE-HK.COM	

TRANSPORT

物流运输



Type 型号	SHZ-95B立式水循环真空泵
Package Size(mm) 尺寸	450*340*840
Weight(Kg) 重量(Kg)	35Kg



PACKAGING:

1. Cartons + wooden cases
2. Custom packaging is available.

打包:

1. 纸箱+木箱。
2. 可以提供定制包装。

SHIPMENT:

1. We can deliver the goods to you by international express, such as DHL, UPS, TNT, EMS and so on.
2. You can also choose the right one according to the schedule and budget.
3. You may choose to use your own shipping agent.

运输:

1. 我们可以通过国际快递将货物运送给您，如DHL，UPS，TNT，EMS等。
2. 您还可以根据时间表和预算选择一个合适的。
3. 您可以选择使用自己的运输代理商。

FAQ

常见问题解答

1. Question: How to choose suitable equipment ?

Answer: You tell us your requirements (or example: capacity, temperature range, pressure, etc), we will recommend you the best suitable equipment.

2. Question: What's the minimum order quantity?

Answer: One set at present.

3. Question: What kind of payment term do you accept?

Answer: Normally we accept TT, Paypal, Western Union, Alibaba Escrow, and L/C. The exactly term with be discussed during the deal.

4. Question: How will you delivery these goods?

Answer: It will be decided according to your quantity. Normally we will choose International Express, Air and sea.

Packaging & Shipping

1.问题：如何选择合适的设备？

答：您告诉我们您的要求（或示例：容量，温度范围，压力等），我们将为您推荐最合适的设备。

2.问题：最低订货量是多少？


答：目前有一套。

3.问题：您接受什么样的付款条件？

答：通常我们接受TT, Paypal, 西联汇款, 阿里巴巴托管和信用证。在交易期间将讨论的确切术语。

4.问题：您将如何交付这些货物？

答：将根据您的数量来决定。通常我们会选择国际快递，空运和海运。
包装和运输

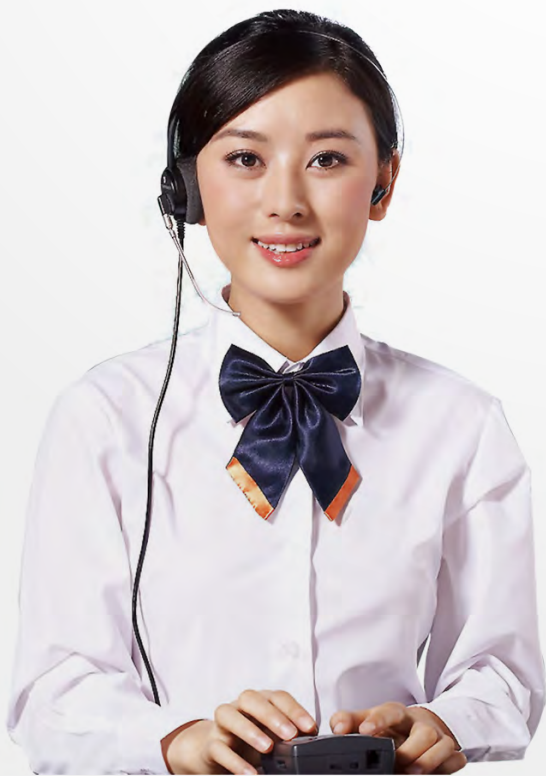

PLEASE CONTACT US

Henan Lanphan Industry Co., Ltd

☎ : +86 - 13783661974

☎ : +86 - 371 - 67447888

☎ : +86 - 371 - 67447999

🌐 : www.rotovap.cn

@ : rotovap@lanphan.com